

Open Access Full Text Article

ORIGINAL RESEARCH

Satisfaction, preference and error occurrence of three dry powder inhalers as assessed by a cohort naïve to inhaler operation

Kwun Nok Mimi Man^{1,*} Zhipeng Tian^{1,*} David Chi-Leung Lam² Jennifer Man Fan Wan³ Kian Cheng Tan-Un¹

¹University of Hong Kong, School of Professional and Continuing Education, Hong Kong SAR, China; ²University of Hong Kong, Department of Medicine, Li Ka Shing Faculty of Medicine, Queen Mary Hospital, Hong Kong SAR, China; ³University of Hong Kong, School of Biological Sciences, Hong Kong

*These authors contributed equally to this work

Background: Inhaled medication is central to the treatment of COPD. Various types of inhaler devices, which directly deliver medication to the lung, have been developed. However, patients often exhibit incorrect techniques of inhaler usage. Effectiveness of therapy may be affected by the ease of device usage, size, convenience of use, durability, clarity of instructions and device preferences of patients. This study compares the satisfaction and preference, as well as error occurrence, with the use of Genuair®, Ellipta™ and Breezhaler™ by healthy subjects in

Subjects and methods: One hundred and thirty healthy Hong Kong Chinese subjects aged ≥40 years without a previous diagnosis of COPD and asthma and with no experience of using dry powder inhalers (DPIs) were recruited. Subjects learned to use the three DPIs by initially reading the instructions and then observing a demonstration with verbal explanation. The number of errors committed was evaluated. Subjects also completed a questionnaire to indicate their satisfaction and preference.

Results: The satisfaction score of comfort for Breezhaler was significantly higher than that for Ellipta ($p \le 0.05$), while the satisfaction score on confidence to have inhaled the entire dose was highest for Genuair compared with Ellipta ($p \le 0.0001$) or Breezhaler ($p \le 0.05$). The overall satisfaction score was significantly higher for Genuair than Ellipta (p≤0.05) or Breezhaler $(p \le 0.01)$. After reading the instructions, the highest number of subjects committing one or more critical errors was with Breezhaler (97) followed by Genuair (70) and then Ellipta (33). Demonstration reduced the number of critical errors made by subjects for each DPI to one third or lower.

Conclusion: Breezhaler seemed to be more comfortable and easy to carry, but users made less critical errors when using Ellipta after reading the instructions only. Genuair provided the clearest indication of correct dose preparation and inhalation.

Keywords: chronic obstructive pulmonary disease, inhalation device, Genuair®, ElliptaTM, BreezhalerTM

Introduction

Inhaled medications, including short-acting and long-acting bronchodilators as well as inhaled corticosteroids, are central to the management of COPD. Various types of inhalers have been developed to enable direct delivery of inhaled medications to the lung, thus minimizing unwanted systemic effects and allowing smaller doses to be used with faster onset of drug action; 1-4 for instance, the pressurized metered-dose inhalers, soft-mist inhalers, dry powder inhalers (DPIs) and nebulizers. DPIs and soft-mist inhalers represent significant improvements from pressurized metered-dose inhalers,

Correspondence: Kian Cheng Tan-Un University of Hong Kong, School of Professional and Continuing Education, 13/F, Fortress Tower, 250 King's Road, North Point, Hong Kong SAR, China Tel +852 2 975 5881 Fax +852 2 707 3620 Email kctan@hkuspace.hku.hk

as they do not contain propellant gases and importantly are breath-actuated, circumventing the need for coordination of inhaler actuation and inspiration by the patient.^{3,5,6} However, only 40%–60% of COPD patients adhered to their prescribed therapy⁷ and, depending on the types of inhalers, between 4% and 94% of patients exhibited incorrect techniques of inhaler device usage.8 Although patients claimed that their techniques of inhaler use is adequate, as many as 94% still committed at least one error when they were asked to demonstrate the correct techniques,9 showing the need for patient education and critical assessments on inhaler device usage by health care professionals. Moreover, the risk of COPD increases by 5-fold for people aged over 65 years as compared to those under 40 years of age. 10 Elderly patients are more likely to have reduced dexterity and cognitive function, which limit their ability to correctly operate inhalers. Inadequate adherence to prescribed regimens and failure to use the correct technique result in reduced effectiveness of therapy. 11,12 Factors affecting patient satisfaction and preference of inhalers include ease of use, size, convenience, durability and clarity of instructions. 13,14

The three most recently launched DPIs for COPD were studied – Genuair® (AstraZeneca, Cambridge, UK), ElliptaTM (GlaxoSmithKline, Brentford, UK) and BreezhalerTM (Novartis, Basel, Switzerland). Genuair is a multidose, mediumresistance DPI for the delivery of aclidinium bromide, a longacting muscarinic agonist (LAMA), or aclidinium bromide in combination with formoterol fumarate dihydrate, a longacting beta agonist (LABA). It provides both visual and acoustic feedback to patients to indicate correct inhalation of the dose. Ellipta is also a multidose DPI for the delivery of umeclidinium bromide (LAMA) with vilanterol (LABA), or fluticasone furoate, an inhaled corticosteroid, and vilanterol

(LABA), or umeclidinium bromide (LAMA) alone. ¹⁷ Breezhaler is a single-dose capsule inhaler for the delivery of glycopyrronium bromide (LAMA) with indacaterol (LABA), or glycopyrronium bromide (LAMA) alone, or indacaterol (LABA) alone. Inhalation causes the capsule to vibrate in the Breezhaler device and such vibrations act as a form of acoustic feedback for the release of medication from the device.

The objectives of this study were to compare the satisfaction and preference, as well as error occurrence with use of these three DPI devices by healthy subjects in Hong Kong.

Subjects and methods Subjects and study design

This was a single-centred, randomized, closed-labeled study. One hundred and thirty healthy Hong Kong Chinese subjects aged ≥40 years who had no previous diagnosis of COPD or asthma and no experience of using DPIs were recruited to the study, thus ruling out any bias on inhaler operation. Subjects who had any condition which impaired their ability to operate the inhaler, who had any history of chronic respiratory diseases, and who were unable to read the instructions, answer subject-reported questionnaires or give written informed consent were excluded from the study. Subjects were randomized with respect to the sequence of learning and the use of the three placebo DPI devices, namely, Genuair, Ellipta and Breezhaler, such that each subject undertook one of the sequences depicted in Figure 1.

Placebo devices of Genuair and Ellipta did not contain any active medication or powder, and in the case of the Breezhaler, empty capsules were used in place of medication-containing capsules. The brand names of the inhalers were not disclosed to the subjects. Subjects were first asked to read a set of written instructions for one of the DPIs and then demonstrate the

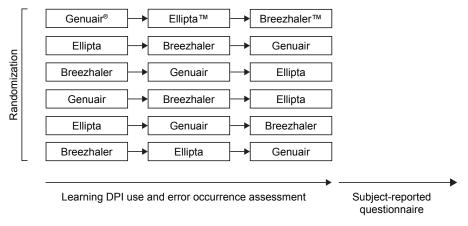


Figure I A schematic diagram showing randomization and sequence of learning to use the different DPI devices. Abbreviation: DPI, dry powder inhaler.

use of the DPI to an assessor. The demonstration was evaluated by the assessor who noted the errors committed based on a checklist (Table S1). Subjects were then given a physical demonstration by the assessor on how to use the DPI with a pre-recorded verbal explanation. Subjects were asked to demonstrate the use of the DPI for a second time, while the assessor observed and noted any error committed. The two learning stages were repeated for the remaining two inhalers. The errors were classified into critical and non-critical errors (Table S1). Critical errors were defined as errors which prevented the subjects from inhaling any of the medications contained in the DPI. Non-critical errors were defined as errors which caused suboptimal administration of the drug or use of the DPI, although the subject could still inhale part of the dose. One main assessor and two assistants were employed to evaluate the errors. The main assessor demonstrated the borderline actions between the correct actions and errors to the assistant. All of them assessed several subjects together at the beginning of the project to ensure the standard was agreed. Then the three assessors evaluated the errors individually. After the subjects had demonstrated the use of all three DPIs, they completed a questionnaire to indicate their level of satisfaction and preference for various attributes of the three inhalers, using a Likert-type scale with a score from 1 to 7. Subjects were also asked, using a 100-point scale for overall satisfaction, to indicate how willing they would be to use each DPI if they were diagnosed with COPD.

The study protocol was approved by the Research Ethics Subcommittee of the University of Hong Kong School of Professional and Continuing Education. Written consent was provided by each subject prior to the commencement of the study procedure.

Statistical analyses

Age of subjects and satisfactory scores were presented as mean \pm standard error of the mean. The occurrence of one or more critical errors during the use of each DPI was analyzed by McNemar's test. Briefly, the three DPIs were compared in a pair-wise manner, that is, Genuair vs Ellipta, Ellipta vs Breezhaler and Genuair vs Breezhaler. In each comparison between two DPIs A vs B, subjects were classified into four groups: 1) committing critical errors when using both A and B, 2) committing critical errors when using A only (the number of which is designated a), 3) committing critical errors when using B only (designated b) and 4) committing no errors when using A or B. The chi-square statistics was calculated by the formula $(a-b)^2/(a+b)$.

Satisfaction scores for each attribute and overall satisfaction scores for each DPI were analyzed for statistical

significance by one-way analysis of variance followed by post hoc Tukey's test. For comparisons of scores for each attribute between male and female subjects and between subjects < 60 and ≥60 years of age, two-way analysis of variance followed by Sidak's multiple comparisons test was used. Minimally important difference was analyzed by comparing the differences in satisfactory score for each attribute between DPIs with the SDs of the scores, with the differences >0.2 SD and <0.5 SD being a small effect and \ge 0.5 SD and <0.79 SD being a medium effect. 18 Distribution of DPI preference was analyzed by the procedure described in Sharpe's study¹⁹ for deviation from uniform distribution by the chi-square goodness-of-fit tests. In order to determine the DPI causing the deviation, we computed the standardized residuals, R. When the absolute value of the standardized residual exceeds 2, this shows lack of fit in the cells, with $p \le 0.05$.²⁰

Results

Subjects

Data on satisfaction, preference and error occurrence for the three DPIs were collected from a total of 130 healthy Chinese subjects (aged 41–84 years, mean age=59.82±0.88 years) consisting of 61 males (aged 42–84 years, mean age=60.98±0.93 years) and 69 females (aged 41–79 years, mean age=58.78±0.82 years).

Satisfaction

Subjects were asked to evaluate and indicate their satisfaction on various attributes of the three DPIs using a seven-point Likert-type scale (Figure 2). In terms of comfort of the DPIs, the satisfaction score for Breezhaler was significantly higher than that for Ellipta $(5.13\pm0.13 \text{ vs } 4.70\pm1.65, p=0.047)$. However, the satisfaction score for Breezhaler in terms of ease of dose preparation was lowest among the DPIs, being significantly lower than those of both Ellipta $(4.75\pm0.15 \text{ vs } 5.24\pm0.15, p=0.0463)$ and Genuair $(4.75\pm0.15 \text{ vs } 5.39\pm0.13, p=0.0038)$.

When being asked about the clarity of the DPIs to indicate correct dose preparation, subjects indicated significantly higher satisfaction for Genuair than Ellipta $(5.69\pm0.12 \text{ vs } 5.25\pm0.15, p=0.0120)$. In terms of clarity of DPIs to indicate correct inhalation, subjects had greater satisfaction with Genuair compared with Ellipta $(5.65\pm0.13 \text{ vs } 5.11\pm0.15, p=0.0014)$. The satisfaction score on confidence to have inhaled the entire dose when using the DPIs was highest for Genuair compared with Ellipta $(5.43\pm0.13 \text{ vs } 4.62\pm0.15, p\leq0.0001)$ and Breezhaler $(5.43\pm0.13 \text{ vs } 4.98\pm0.13, p=0.0328)$.

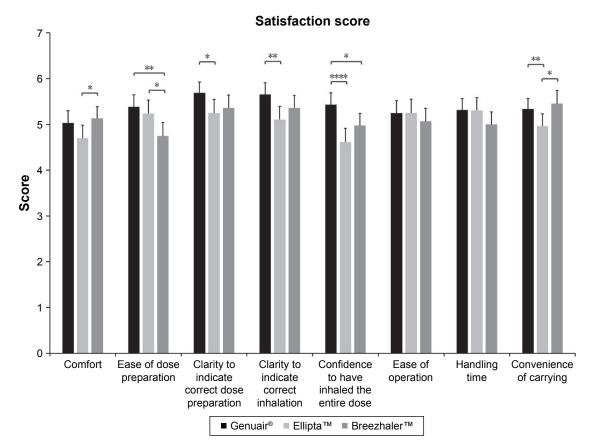


Figure 2 Satisfaction scores on different attributes of the three DPIs using a seven-point Likert-type scale. **Notes:** Mean values are plotted and the error bars represent the SEM. * $p \le 0.05$; *** $p \le 0.01$; ***** $p \le 0.0001$. **Abbreviations:** DPIs, dry powder inhalers; SEM, standard error of the mean.

No significant difference was observed in the satisfaction score for ease of operation and handling time among the three DPIs (Figure 2). Finally, subjects considered Breezhaler $(5.45\pm0.15 \text{ vs } 4.96\pm0.14, p=0.0304)$ and Genuair $(5.34\pm0.12 \text{ vs } 4.96\pm0.14, p=0.0087)$ to be more convenient to carry than Ellipta. Similar trends were seen in the subanalysis by gender and age (younger than 60 years of age vs 60 years of age and above) and these data can be found in Figures S1 and S2. An analysis that directly compares the effects of gender (male vs female) and age ($<60 \text{ vs } \ge 60 \text{ years of age}$) on the satisfactory scores of the different attributes is provided in Figures S3 and S4, respectively.

The satisfaction score data were further analyzed to determine whether the effect sizes exceed the threshold for clinical relevance by considering the minimally important difference using a distribution-based method. Conventional benchmarks were used, wherein an effect size > 0.2 SD and < 0.5 SD is a small effect and an effect size ≥ 0.5 SD and < 0.79 SD is a medium effect. After analyzing the differences between the satisfaction scores of each DPI for each attribute (Tables 1 and S2), we found that all effects which reached statistical significance were small effects, with

the exception of clarity to indicate correct inhalation, with Genuair being more superior than Ellipta and approaching a medium effect.

The overall satisfaction score, used as a measure to indicate the willingness of subjects to use each of the DPIs

Table I Differences of satisfactory scores between each DPI for each attribute

Attribute	Genuair®-	Ellipta-	Genuair-
	Ellipta™	Breezhaler™	Breezhaler
Comfort	0.3308	0.4308 ^a	0.1000
Ease of dose preparation	0.1462	0.4923 ^a	0.6385ª
Clarity to indicate	0.4462ª	0.1154	0.3308
correct dose preparation			
Clarity to indicate	0.5462a	0.2538	0.2923
correct inhalation			
Confidence to have	0.8154 ^b	0.3615	0.4538^{a}
inhaled the entire dose			
Ease of operation	0.0077	0.1846	0.1769
Handling time	0.0077	0.3077	0.3154
Convenience of carrying	0.3769^{a}	0.4923ª	0.1154
Overall satisfaction score	5.5615ª	2.2615	7.823 la

Note: ³Indicates significance (see Figure 2) and a small effect, while ⁵Indicates significance and an effect approaching medium.

Abbreviation: DPI, dry powder inhaler.



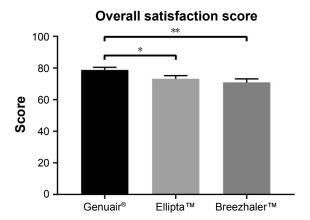


Figure 3 Overall satisfaction scores on a 0–100 scale to indicate the willingness of subjects to use each DPI (*p=0.05; **p=0.01). **Abbreviation:** DPI, dry powder inhaler.

if they were diagnosed with COPD, was significantly higher for Genuair than for both Ellipta (78.79 \pm 1.70 vs 73.23 \pm 1.94, p=0.0132) and Breezhaler (78.79 \pm 1.70 vs 70.97 \pm 2.15, p=0.0063; Figure 3).

Preferences

Subjects were asked to indicate their preferred inhaler for each attribute evaluated (Table 2). The distribution of preference for each attribute was analyzed for significant deviations from uniform distribution (ie, 33.3% for each DPI) using chi-square tests.¹⁷

Genuair was the preferred DPI in terms of clarity to indicate correct dose preparation (38.5%, R=2.137), clarity to indicate correct inhalation (37.7%, R=2.299) and confidence to have inhaled the entire dose (42.3%, R=2.824). In terms of confidence to have inhaled the entire dose, the standardized residual for Ellipta was -2.716, indicating a significant lack of preference (16.2%). Breezhaler was the preferred DPI for convenience of carrying, as preferred by 42.3% of subjects (R=2.372), while Ellipta was less preferred by subjects for this attribute (16.2%, R=-3.004). For the attributes of

comfort, ease of dose preparation, ease of operation and handling time, the distributions were not found to deviate significantly.

Error occurrence during use of DPIs

After reading the instructions, the highest number of subjects committing one or more critical errors was with Breezhaler (96 [73.8%]) followed by Genuair (70 [53.8%]) and then Ellipta (33 [25.4%]), as shown in Figure 4A. The likelihood for subjects to commit critical errors when using Breezhaler was higher than those when using either Genuair ($p \le 0.0001$, chi-square statistic=15.36) or Ellipta ($p \le 0.0001$, chi-square statistic=55.90). Similarly, subjects were more likely to commit one or more critical errors when using Genuair than Ellipta ($p \le 0.0001$, chi-square statistic=27.94; Figure 4B).

For all three DPIs, the number of subjects who committed critical errors reduced after visual demonstration with verbal instructions. The highest number of subjects committing one or more critical errors, after the demonstration, was with Breezhaler (30 [23.1%]) followed by Ellipta (12 [9.2%]) and then Genuair (11 [8.5%]), as shown in Figure 4A. The likelihood for subjects to commit critical errors when using Breezhaler was still higher than for those using either Genuair ($p \le 0.001$, chi-square statistic=12.45) or Ellipta ($p \le 0.01$, chi-square statistic=10.80). There was no significant difference between Genuair and Ellipta (p = 0.8084, chi-square statistic=0.0588; Figure 4B).

The number and percentage of subjects who committed each non-critical error is depicted in Table 3. The total number of non-critical errors committed during use of each DPI is presented in Table S3.

Discussion


The correct and preferred use of DPIs is essential to deliver the relevant drugs appropriately to the patient. Besides the efficacy of bronchodilators and corticosteroids,

Table 2 Preference of subjects for different attributes of the three DPIs

Attributes	p-value of chi-	% Preference (st	andardized residual,	R)	
	square goodness- of-fit test	Genuair [®]	Ellipta™	Breezhaler™	
Comfort	0.097	33.8 (0.801)	21.5 (-1.761)	34.6 (0.971)	
Ease of dose preparation	0.569	29.2 (-0.265)	34.6 (0.847)	27.7 (-0.582)	
Clarity to indicate correct dose preparation ^a	0.014	38.5 (2.137)	19.2 (-1.973)	27.7 (-0.164)	
Clarity to indicate correct inhalation ^a	0.013	37.7 (2.299)	19.2 (-1.738)	24.6 (-0.561)	
Confidence to have inhaled the entire dose ^a	<0.0001	42.3 (2.824)	16.2 (-2.716)	28.5 (-0.109)	
Ease of operation	0.294	33.1 (0.474)	34.6 (0.791)	24.6 (-1.265)	
Handling time	0.116	32.3 (0.764)	33.1 (0.927)	20.8 (-1.691)	
Convenience of carrying ^a	0.001	33.8 (0.632)	16.2 (-3.004)	42.3 (2.372)	

Note: *Significant deviation from uniform distribution, with $p \le 0.05$, R > 2 or R < -2.

Abbreviation: DPIs, dry powder inhalers.

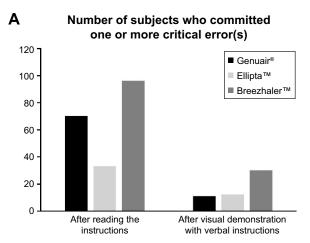

B Number of subjects who committed one or more critical error(s) during use of one of the three DPIs after reading the instructions

		Elli	pta			Bree	zhaler			Breez	haler
		Yes	No			Yes	No			Yes	No
Genuair	Yes	27	43	Ellipta	Yes	29	4	Genuair	Yes	61	9
	No	6	54		No	67	30		No	35	25
p-value:			1.25e-07*	p-value:			7.62e-14*	p-value:			8.87e-5*
Chi-square	e statistic:		27.939	Chi-square	e statistic:		55.901	Chi-square	e statistic:		15.364

After visual demonstration with verbal instructions

		Ellip	ota			Breez	haler			Breez	haler
		Yes	No			Yes	No			Yes	No
Genuair	Yes	3	8	Ellipta	Yes	6	6	Genuair	Yes	6	5
	No	9	110		No	24	94		No	24	95
p-value:			0.8084	p-value:			0.00102*	p-value:			0.00042*
Chi-square	statistic:		0.0588	Chi-square	statistic:		10.8	Chi-square	statistic:		12.448

Figure 4 (A) Number of subjects who committed one or more critical errors during the two learning stages. (B) McNemar's tests for the number of subjects who committed one or more critical errors during the two learning stages.

Note: **p*≤0.05

Abbreviation: DPIs, dry powder inhalers.

patient preference and the skills required in using inhaler devices can affect the effectiveness of the treatment. It is important that an inhaler device is easy to use, so that the correct dosage is administered. In this study, 130 healthy subjects were recruited and the numbers of male (n=61) and female (n=69) subjects were evenly distributed. The subjects were aged 41-84 years with the average age 59.82±0.88 years. As the prevalence of COPD is higher in people aged 60 years or above compared to those aged 40-60 years, 22,23 the age distribution of our subjects also covered a cross section of the high and low prevalence groups. Besides disease conditions such as Huntington's chorea, stroke and so on, which cause dysfunction of brain, 24,25 degeneration of brain structure with aging is highly related to reduced dexterity and cognitive function. 26,27 Thus, ease of use is very important for elderly patients. The study did not include the education level of the subjects. Although education background might affect the ability to understand instructions, Blasi et al reported no difference in the time required to perform a correct inhalation of Genuair between different education levels.²⁸

In our study, subjects rated Breezhaler with the highest score for "comfort" (Figure 2). Interestingly, the score given by males did not show any significant difference among the three DPIs, and the high score of Breezhaler for "comfort" was mainly contributed by females (Figure S2). This could be related to the shape and size of Breezhaler. Compared to the other two DPIs, Breezhaler is smaller and the medicine is not built into the device, which could also explain the higher score for "convenience of carrying". However, the design seems to also affect dose preparation, as the user required extra steps to put a capsule into Breezhaler. The subjects in our study gave a significantly lower score for Breezhaler for "ease of dose preparation".

Both Ellipta and Genuair are equipped with acoustic feedback mechanism to indicate that a dose has been correctly prepared. Genuair also has a visual feedback mechanism in the form of a red/green window to indicate empty and loaded dose. These visual and acoustic feedback mechanisms also indicate correct inhalation of a dose when enough inhalation flow is generated in Genuair. The higher scores for Genuair in "clarity to indicate correct dose preparation" and

 Table 3 Number and percentage of subjects who committed non-critical errors

Non-critical errors	After	After visual
	reading the	demonstration with
	instructions	verbal instructions
Genuair®		
Did not hold the inhaler at the correct position during dose preparation	34 (26.15)	4 (3.08)
Did not breathe out completely before inhalation ^a	52 (40.00)	19 (14.62)
Breathed into the inhaler before inhalation	11 (8.46)	4 (3.08)
Did not hold the inhaler at the correct position during inhalation	2 (1.54)	0 (0.00)
Did not breathe in strongly and deeply during inhalation ^a	52 (40.00)	5 (3.85)
Stopped breathing in when a "click" sound was heard	9 (6.92)	14 (10.77)
Pressed the button during or after inhalation	29 (22.31)	14 (10.77)
Did not hold breath for long enough after inhalation ^a	49 (37.69)	20 (15.38)
Breathed out quickly and strongly after holding the breath/breathed out through the mouth	32 (24.62)	28 (21.54)
Did not replace the protective cap	10 (7.69)	4 (3.08)
Ellipta™		
Opened the cover and closed it without inhalation	5 (3.85)	3 (2.31)
Did not hold the inhaler at the correct position during dose preparation	15 (11.54)	8 (6.15)
Shook the inhaler after loading a dose	9 (6.92)	0 (0.00)
Did not breathe out completely before inhalation ^a	66 (50.77)	21 (16.15)
Breathed into the inhaler before inhalation	9 (6.92)	7 (5.38)
Blocked the air vent with fingers during inhalation	14 (10.77)	2 (1.54)
Did not hold the inhaler at the correct position during inhalation	15 (11.54)	11 (8.46)
Did not breathe in strongly and deeply during inhalation	49 (37.69)	27 (20.77)
Did not hold breath for long enough after inhalation ^a	65 (50.00)	25 (19.23)
Breathed out quickly and strongly after holding the breath ^a	52 (40.00)	23 (17.69)
Breezhaler™		
Did not close the inhaler properly after placing the capsule (with a "click" sound)	I (0.77)	I (0.77)
Did not hold the inhaler at the correct position during dose preparation	7 (5.38)	7 (5.38)
Pressed the two side buttons more than once	33 (25.38)	11 (8.46)
Took inhalation posture without releasing the side buttons fully	18 (13.85)	6 (4.62)
Did not breathe out completely before inhalation ^a	58 (44.62)	17 (13.08)
Breathed into the inhaler before inhalation	4 (3.08)	6 (4.62)
Did not hold the inhaler at the correct position during inhalation	10 (7.69)	7 (5.38)
In case a whirring noise was not heard, did not open the inhaler to loosen the capsule ^a	37 (28.46)	15 (11.54)
When checking a stuck capsule, did not loosen the capsule by tapping the base of the inhaler	4 (3.08)	2 (1.54)
Did not hold breath after inhalation for long enough ^a	42 (32.31)	16 (12.31)
Breathed out quickly and strongly after holding the breath	33 (25.38)	20 (15.38)
Removed the empty capsule by touching the capsule in the inhaler	8 (6.15)	4 (3.08)
Did not replace the cap	17 (13.08)	8 (6.15)

Notes: Data show the number (percentage) of subjects. ^aThe three most common non-critical errors for each DPI are indicated. **Abbreviation:** DPI, dry powder inhaler.

"clarity to indicate correct inhalation" (Figure 2) reflected that subjects were more satisfied with the device's feedback mechanism. The subjects also rated a significantly higher score for Genuair in "confidence to have inhaled the entire dose". Rajan and Gogtay suggested that confidence was important to new users of DPI, as this would lead to better adherence to therapy.²⁹ The clearer indications of dose preparation and correct inhalation together with higher confidence could explain the significantly higher overall satisfaction score for Genuair. A previous study in Germany, Spain and the UK compared the satisfaction of COPD patients on Genuair and Breezhaler.³⁰ The results also showed that a higher overall satisfaction score was found in Genuair than Breezhaler.

The satisfaction scores of different attributes were analyzed by gender and age (Figures S1 and S2). The pattern of the scores for the three DPIs was similar between males and females. When analyzing the satisfaction scores according to age groups, the patterns were also very similar. However, in terms of "ease of operation" and "handling time", the younger group (aged 40–60) reported a significantly higher score for Ellipta than Breezhaler, while the older group (aged 60 and above) gave a higher score to Breezhaler. The satisfaction scores of different attributes were further compared by male vs female and the younger group vs older group (Figures S3 and S4). The only significant result was found in "ease of operation", where the older group gave a significantly lower score to Ellipta. These results suggested

that the handling procedures of Ellipta may be difficult for the older group.

The errors during the use of DPIs were also evaluated based on the list of errors (Table S1). The handling procedures and complexity of the three DPIs are different, and thus, the error list is designed based on the instruction manual of each DPI. After reading the instructions only, subjects were less likely to commit critical errors when using Ellipta than Genuair or Breezhaler (Figure 4B). The data suggest that Ellipta may be easier to learn by reading the instructions alone. Regardless of reading the instructions only or after observation of the visual demonstration with verbal instructions, the number of subjects committing critical errors was higher with Breezhaler. These results may indicate that the design of Breezhaler is more complex, especially for new users. Interestingly, subjects also gave a lower satisfaction score for Breezhaler on "ease of operation" and "handling time" (Figure 2). A previous study in Japan also showed that the volunteers made more errors when using Breezhaler compared to Ellipta.¹⁷ It is also worth noting that visual demonstration with verbal instructions reduced the number of subjects making critical errors for each DPI to one third or lower compared to reading the instructions only, with the effect being the greatest for Genuair. This finding suggests that visual demonstration with verbal instructions is crucial to new users of DPIs. In a study conducted by Svedsater et al, subjects also reported that with briefing and demonstration, they could understand how to use DPIs even without reading the instructions.31

This study excluded subjects who had been diagnosed with COPD or asthma and had experience of using DPIs. Thus, the results would represent newly diagnosed COPD patients who are device naïve. In the study of the usability of Genuair, Blasi et al reported that there was no difference in the time required to perform a correct inhalation when considering the different previous experience of inhaler devices. However, patients' experience of using DPI may affect their satisfaction. The results could be used as reference for future validation studies in COPD patients with previous experience in using DPIs. A dose counter is found in Genuair and Ellipta to indicate the number of doses left in the device. The satisfaction to this function was not tested, as the subjects demonstrated the use of each DPI for only two times.

Conclusion

The results of our study showed that each DPI has its own user advantages: Breezhaler seemed to be more comfortable to use and easy to carry. Device-naive users made less critical

errors when using Ellipta after reading the instructions only. Genuair gave significantly better indications of correct dose preparation and dose inhalation. Overall, healthy volunteers had the greatest satisfaction scores when using Genuair. This is important, as satisfaction might improve patient compliance to inhaler therapy. Instead of learning to use a DPI by reading the instructions only, visual demonstration with verbal instructions given by health care professionals is crucial as they can reduce the critical errors during inhaler use. These results cannot be applied at the individual level, so good practice should be followed.

Acknowledgments

We would like to thank Mr Kwok Yan Ng, University of Hong Kong, School of Biological Sciences for his help in recruiting participants. We appreciate the support from Dr Hiu Fai Law with the statistical analysis of the data.

Disclosure

The authors report no conflicts of interest in this work.

References

- 1. Barnes PJ. Introduction: how can we improve asthma management? Curr Med Res Opin. 2005;21(Suppl 4):S1–S3.
- Cross S. Asthma inhalation delivery systems: the patient's viewpoint. J Aerosol Med. 2001;14(Suppl 1):S3–S7.
- Magnussen H. Novolizer: how does it fit into inhalation therapy? Curr Med Res Opin. 2005;21(Suppl 4):S39–S46; discussion S47.
- 4. Melani AS. Inhalatory therapy training: a priority challenge for the physician. *Acta Biomed*. 2007;78(3):233–245.
- Dal Negro RW, Povero M. Acceptability and preference of three inhalation devices assessed by the Handling Questionnaire in asthma and COPD patients. *Multidiscip Respir Med.* 2015;11:7.
- Virchow JC. Guidelines versus clinical practice-which therapy and which device? Respir Med. 2004;98(Suppl B):S28–S34.
- Restrepo RD, Alvarez MT, Wittnebel LD, et al. Medication adherence issues in patients treated for COPD. *Int J Chron Obstruct Pulmon Dis*. 2008;3(3):371–384.
- Lavorini F, Magnan A, Dubus JC, et al. Effect of incorrect use of dry powder inhalers on management of patients with asthma and COPD. Respir Med. 2008;102(4):593–604.
- Souza ML, Meneghini AC, Ferraz E, Vianna EO, Borges MC. [Knowledge of and technique for using inhalation devices among asthma patients and COPD patients]. *J Bras Pneumol*. 2009;35(9):824–831. Portuguese.
- Raherison C, Girodet PO. Epidemiology of COPD. Eur Respir Rev. 2009; 18(114):213–221.
- Chrystyn H, Small M, Milligan G, Higgins V, Gil EG, Estruch J. Impact of patients' satisfaction with their inhalers on treatment compliance and health status in COPD. *Respir Med*. 2014;108(2):358–365.
- 12. Lareau SC, Yawn BP. Improving adherence with inhaler therapy in COPD. Int J Chron Obstruct Pulmon Dis. 2010;5:401–406.
- Hodder R, Price D. Patient preferences for inhaler devices in chronic obstructive pulmonary disease: experience with Respimat Soft Mist inhaler. *Int J Chron Obstruct Pulmon Dis*. 2009;4:381–390.
- van der Palen J, Ginko T, Kroker A, et al. Preference, satisfaction and errors with two dry powder inhalers in patients with COPD. Expert Opin Drug Deliv. 2013;10(8):1023–1031.

- Chrystyn H, Niederlaender C. The Genuair[®] inhaler: a novel, multidose dry powder inhaler. *Int J Clin Pract*. 2012;66(3):309–317.
- Magnussen H, Watz H, Zimmermann I, et al. Peak inspiratory flow through the Genuair inhaler in patients with moderate or severe COPD. Respir Med. 2009;103(12):1832–1837.
- Komase Y, Asako A, Kobayashi A, Sharma R. Ease-of-use preference for the ELLIPTA® dry powder inhaler over a commonly used singledose capsule dry powder inhaler by inhalation device-naïve Japanese volunteers aged 40 years or older. *Int J Chron Obstruct Pulmon Dis*. 2014;9:1365–1375.
- Kozma CM, Slaton TL, Monz BU, Hodder R, Reese PR. Development and validation of a patient satisfaction and preference questionnaire for inhalation devices. *Treat Respir Med.* 2005;4(1):41–52.
- Sharpe D. Your chi-square test is statistically significant: now what? Pract Assess Res Eval. 2015;20(8):1–10.
- Agresti A. An Introduction to Categorical Data Analysis. Hoboken, New Jersey, USA. Wiley; 2007.
- Guyatt GH, Osoba D, Wu AW, Wyrwich KW, Norman GR; Clinical Significance Consensus Meeting Group. Methods to explain the clinical significance of health status measures. *Mayo Clin Proc.* 2002;77(4): 371–383.
- Menezes AM, Perez-Padilla R, Jardim JR, et al; PLATINO Team. Chronic obstructive pulmonary disease in five Latin American cities (the PLATINO study): a prevalence study. *Lancet*. 2005;366(9500): 1875–1881.
- Global Strategy for Diagnosis, Management, and Prevention of COPD 2016. Global Initiative for Chronic Obstructive Lung Disease – GOLD. Available from: http://goldcopd.org/global-strategy-diagnosis-management-prevention-copd-2016/. Accessed May 26, 2017.

- Montoya A, Price BH, Menear M, Lepage M. Brain imaging and cognitive dysfunctions in Huntington's disease. *J Psychiatry Neurosci*. 2006;31(1):21–29.
- Sunderland A, Bowers MP, Sluman SM, Wilcock DJ, Ardron ME. Impaired dexterity of the ipsilateral hand after stroke and the relationship to cognitive deficit. Stroke. 1999;30(5):949–955.
- Park DC, Reuter-Lorenz P. The adaptive brain: aging and neurocognitive scaffolding. *Annu Rev Psychol*. 2009;60:173–196.
- Seidler RD, Bernard JA, Burutolu TB, et al. Motor control and aging: links to age-related brain structural, functional, and biochemical effects. *Neurosci Biobehav Rev.* 2010;34(5):721–733.
- 28. Blasi F, Canonica GW, Centanni S, et al. Genuair® usability test: results of a national public survey of the elderly. *COPD*. 2016;13(3): 367–371
- Rajan SK, Gogtay JA. Ease-of-use, preference, confidence, and satisfaction with Revolizer(*), a novel dry powder inhaler, in an Indian population. *Lung India*. 2014;31(4):366–374.
- Pascual S, Feimer J, De Soyza A, et al. Preference, satisfaction and critical errors with Genuair and Breezhaler inhalers in patients with COPD: a randomised, cross-over, multicentre study. NPJ Prim Care Respir Med. 2015;25:15018.
- Svedsater H, Dale P, Garrill K, Walker R, Woepse MW. Qualitative assessment of attributes and ease of use of the ELLIPTATM dry powder inhaler for delivery of maintenance therapy for asthma and COPD. BMC Pulm Med. 2013;13:72.

Supplementary materials

 Table ST List of critical and non-critical errors for each dry powder inhaler

Table S1 List of critical and non-critical errors for each dry powder inhaler						
	Errors	Critical (C) or non-critical (NC)				
Genuair®						
1	Did not remove cap before inhalation	С				
2	Did not hold the inhaler at the correct position during dose preparation	NC				
3	Did not press the button all the way down before inhalation	С				
4	Did not release the button before inhalation	С				
5	Inhaled when the control window was still red	С				
6	Did not breathe out completely before inhalation	NC				
7	Breathed into the inhaler before inhalation	NC				
8	Did not close lips firmly around the mouthpiece	С				
9	Did not hold the inhaler at the correct position during inhalation	NC				
10	Did not breathe in strongly and deeply during inhalation	NC				
11	Stopped breathing in when a "click" sound was heard	NC				
12	Pressed the button during or after inhalation	NC				
13	Did not hold breath for long enough after inhalation	NC				
14	Breathed out quickly and strongly after holding breath/breathed out through the mouth	NC				
15	Did not repeat the inhalation even when the control window was green	С				
16	Did not replace the protective cap	NC				
Ellipta™						
1	Could not open the cover without instructions	С				
2	Opened the cover and closed it without inhalation	NC				
3	Did not hold the inhaler at the correct position during dose preparation	NC				
4	Shook the inhaler after loading a dose	NC				
5	No "click" sound after sliding the cover open	С				
6	Did not breathe out completely before inhalation	NC				
7	Breathed into the inhaler before inhalation	NC				
8	Did not close lips firmly around the mouthpiece	С				
9	Blocked the air vent with fingers during inhalation	NC				
10	Did not hold the inhaler at the correct position during inhalation	NC				
П	Did not breathe in strongly and deeply during inhalation	NC				
12	Did not hold breath for long enough after inhalation	NC				
13	Breathed out quickly and strongly after holding breath	NC				
14	Did not close the cover as far as possible	C				
Breezhaler™						
1	Did not remove cap before inhalation	С				
2	Placed a capsule directly into the mouthpiece	C				
3	Did not close the inhaler properly after placing the capsule (with a "click" sound)	NC				
4	Did not press the two side buttons simultaneously until a "click" was heard before inhalation	C				
5	Did not hold the inhaler at the correct position during dose preparation	NC				
6	Pressed the two side buttons more than once	NC				
7	Took inhalation posture without releasing the side buttons fully	NC				
8	Did not breathe out completely before inhalation	NC				
9	Breathed into the inhaler before inhalation	NC				
10	Did not close lips firmly around the mouthpiece	C				
11	Did not close hips in may a build the moduliplece Did not hold the inhaler at the correct position during inhalation	NC				
12	Did not breathe in strongly and deeply during inhalation (no whirring noise was heard)	C				
13	In case a whirring noise was not heard, did not open the inhaler to loosen the capsule	NC				
14	When checking a stuck capsule, did not loosen the capsule by tapping the base of the inhaler	NC				
15	Did not hold breath after inhalation for long enough	NC				
16	Breathed out quickly and strongly after holding breath	NC				
17	Did not remove the empty capsule	C				
18		NC				
19	Removed the empty capsule by touching the capsule in the inhaler Did not replace the cap	NC				

Satisfaction scores according to gender

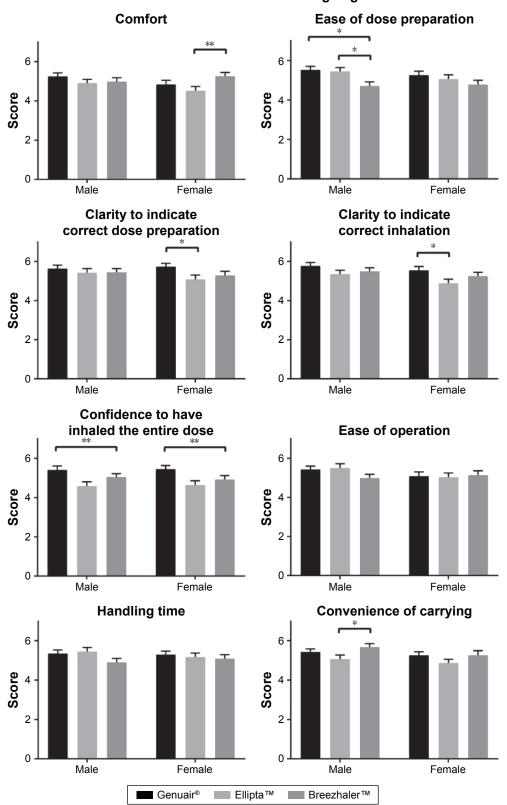


Figure S1 Satisfaction scores on different attributes of the three DPIs as evaluated by male and female subjects. **Notes:** Mean values are plotted and the error bars represent the SEM. * $p \le 0.05$; ** $p \le 0.01$. Abbreviations: DPIs, dry powder inhalers; SEM, standard error of the mean.

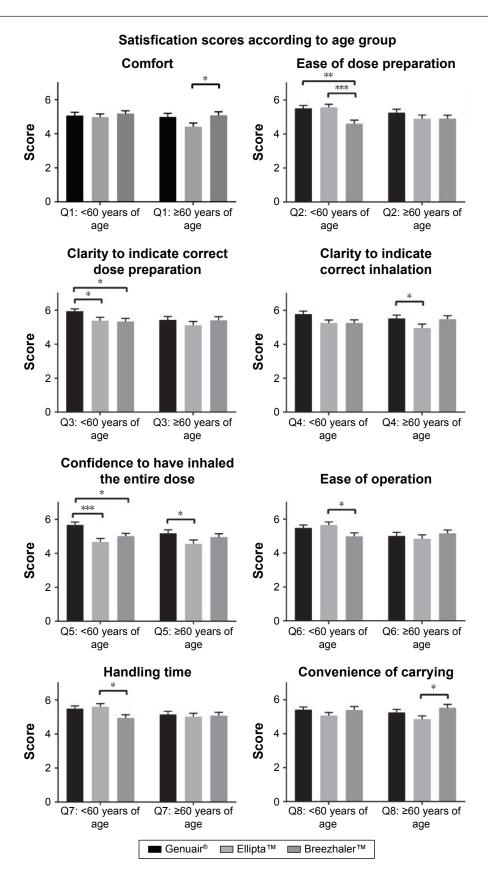


Figure S2 Satisfaction scores on different attributes of the three DPIs as evaluated by subjects below 60 years of age or above or equal to 60 years of age. Notes: Mean values are plotted and the error bars represent the SEM. * $p \le 0.05$; ** $p \le 0.01$; *** $p \le 0.001$. Abbreviations: DPIs, dry powder inhalers; SEM, standard error of the mean.

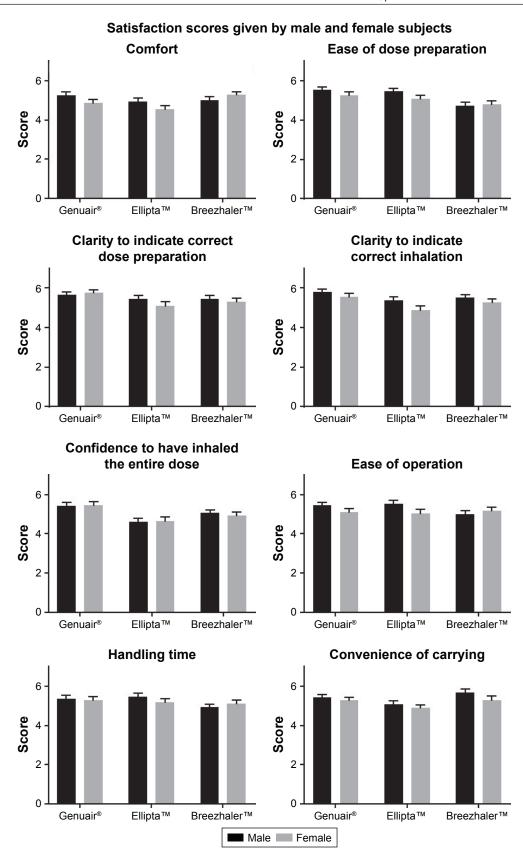
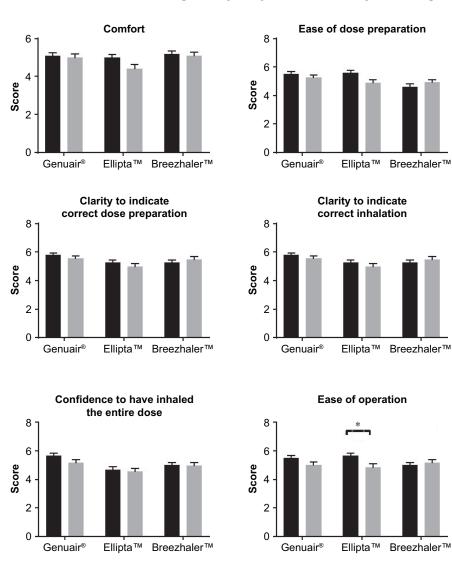


Figure S3 Satisfaction scores on different attributes of the three DPIs by male and female subjects. Notes: Mean values are plotted and the error bars represent the SEM. No significance was detected. Abbreviations: DPIs, dry powder inhalers; SEM, standard error of the mean.

Satisfaction scores given by subjects <60 or ≥60 years of age



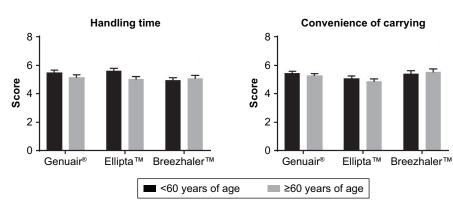


Figure S4 Satisfaction scores on different attributes of the three DPIs by subjects <60 and ≥60 years of age. **Notes:** Mean values are plotted and the error bars represent the SEM. * $p\le0.05$. **Abbreviations:** DPIs, dry powder inhalers; SEM, standard error of the mean.

Table S2 0.2 SD and 0.5 SD of satisfactory scores of the three DPIs

Attribute	Genuair [®]		Ellipta™		Breezhaler™	
	0.2 SD	0.5 SD	0.2 SD	0.5 SD	0.2 SD	0.5 SD
Comfort	0.315	0.787	0.331	0.827	0.302	0.756
Ease of dose preparation	0.304	0.761	0.338	0.844	0.340	0.850
Clarity to indicate correct dose preparation	0.275	0.687	0.349	0.872	0.325	0.812
Clarity to indicate correct inhalation	0.297	0.744	0.331	0.827	0.314	0.785
Confidence to have inhaled the entire dose	0.303	0.757	0.349	0.872	0.305	0.764
Ease of operation	0.317	0.793	0.347	0.868	0.326	0.815
Handling time	0.289	0.723	0.319	0.797	0.315	0.787
Convenience of carrying	0.265	0.663	0.310	0.776	0.338	0.845
Overall satisfaction score	3.870	9.674	4.434	11.084	4.893	12.231

Abbreviation: DPIs, dry powder inhalers.

Table S3 Total number of non-critical errors committed during use of the three DPIs

DPIs	After reading	After visual demonstration
	the instructions	with verbal instructions
Genuair®	280	112
Ellipta™	299	127
Breezhaler™	272	120

Abbreviation: DPIs, dry powder inhalers.

International Journal of COPD

Publish your work in this journal

The International Journal of COPD is an international, peer-reviewed journal of therapeutics and pharmacology focusing on concise rapid reporting of clinical studies and reviews in COPD. Special focus is given to the pathophysiological processes underlying the disease, intervention programs, patient focused education, and self management protocols.

This journal is indexed on PubMed Central, MedLine and CAS. The manuscript management system is completely online and includes a very quick and fair peer-review system, which is all easy to use. Visit http://www.dovepress.com/testimonials.php to read real quotes from published authors.

 $\textbf{Submit your manuscript here:} \ \texttt{http://www.dovepress.com/international-journal-of-chronic-obstructive-pulmonary-disease-journal-obstructive-pulmonary-disease-journal-obstructive-pulmonar$

Dovepress